

SOMNI

SCIENTIFIC



INDUCTION CHAMBERS ACTIVE VACUUM AND PASSIVE STYLES

Induction Chambers

Active Vacuum and Passive Styles

The SOMNI line of induction chambers are designed with slide top lids that allow for one handed operation. Unlike a “hinge lid” the slide top keeps gases from being “drawn” into the workplace when the chamber is opened.

► SPECIFICATIONS

ACTIVE LARGE: 5.10 lbs / 11.5 in L 6 in W 6 in H

ACTIVE MEDIUM: 3.10 lbs / 9.25 in L 5.25 in W 5.25 in H

PASSIVE MEDIUM: 3.10 lbs / 9.25 in L 5.25 in W 5.25 in H

PASSIVE LARGE: 5.10 lbs / 11.5 in L 6 in W 6 in H



QUALITY:



- Constructed of durable, shatter resistant polycarbonate.
- Can be safely cleaned with common clinical disinfectants.

SAFETY:

- Prevents unnecessary exposure to personnel and environment (NIOSH 2ppm) when used with active waste gas removal system EPS-3 or W.A.G. filter (active chamber only).

CONVENIENT:

- One-handed operation
- Color coded inlet and outlet ports.
 - Compatible with 19mm tubing & 15mm adapters
- MRI compatible

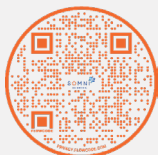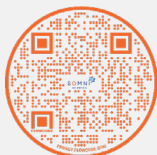

1-877-637-3625

info@somniscientific.com

www.somniscientific.com

Proudly machined, assembled, and calibrated in the USA