

# SOMNI

SCIENTIFIC



## INDUCTION CHAMBERS ACTIVE VACUUM AND PASSIVE STYLES

### Induction Chambers

#### Active Vacuum and Passive Styles

The SOMNI line of induction chambers are designed with slide top lids that allow for one handed operation. Unlike a “hinge lid” the slide top keeps gases from being “drawn” into the workplace when the chamber is opened.

### ► SPECIFICATIONS

**ACTIVE LARGE:** 5.10 lbs / 11.5 in L 6 in W 6 in H

**ACTIVE MEDIUM:** 3.10 lbs / 9.25 in L 5.25 in W 5.25 in H

**PASSIVE MEDIUM:** 3.10 lbs / 9.25 in L 5.25 in W 5.25 in H

**PASSIVE LARGE:** 5.10 lbs / 11.5 in L 6 in W 6 in H



### QUALITY:

- Constructed of durable, shatter resistant polycarbonate.
- Can be safely cleaned with common clinical disinfectants.




### SAFETY:

- Prevents unnecessary exposure to personnel and environment (NIOSH 2ppm) when used with active waste gas removal system EPS-3 or W.A.G. filter (active chamber only).



### CONVENIENT:

- One-handed operation
- Color coded inlet and outlet ports.
  - Compatible with 19mm tubing & 15mm adapters
- MRI compatible

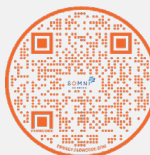

1-877-637-3625

info@somniscientific.com

www.somniscientific.com

Proudly machined, assembled, and calibrated in the USA